MADE IN ITALY

MAGIC CARE®

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MrBeetle

*Congratulations on your purchase and thank you for choosing our unit.
Our goal is to fully satisfy customers by offering them cutting-edge systems
for treating ailments of the upper airways.*



***Carefully read these instructions and keep them in a safe place for future reference.
Use the unit only as described in this instruction manual. This is a medical device for home
use and may only be used with medication prescribed or recommended by a doctor.***

Visit our Internet site www.flaem.it to view the whole range of Flaem products.

G
B

THE UNIT STANDARD FITTINGS INCLUDE:

- (1) Aerosoltherapy unit
- (2) Nebuliser holder
- (3) Air inlet
- (4) Air filter Code. 12960
- (5) On/off switch
- (6) Socket for power cord
- (7) Power cord
- (8) Connection tube
- (9) **RF6 Basic²** Nebuliser (see also page 6)
- (10) Mouthpiece
- (11) Adult mask (PVC)
- (12) Child mask (PVC)
- (13) Carrying bag

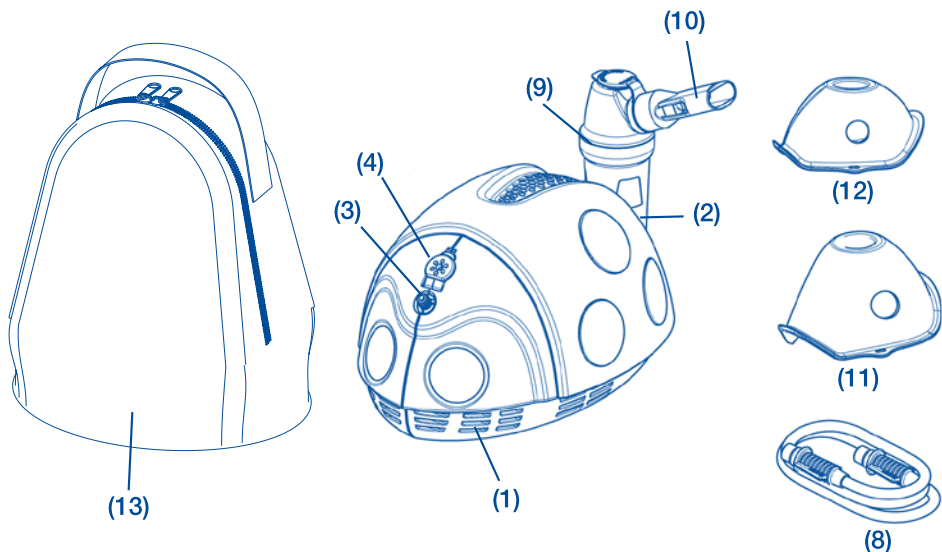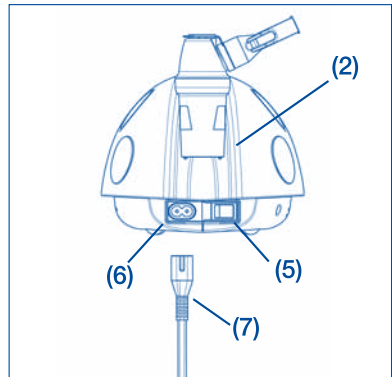TREATMENT OF THE WHOLE RESPIRATORY TRACT

1. INSTRUCTIONS FOR RF6 BASIC² NEBULISER

Prior to use, the nebuliser and the accessories must be cleaned according to the instructions provided in the “Cleaning, disinfection and sterilisation” paragraph.

⚠ Remove the unit from the carrying bag (13) before use!

- 1 Plug the power cord (7) into the socket (6) at the rear of the unit and plug the power cord into a power socket corresponding to the voltage of the unit.
- 2 Carefully wash your hands before going ahead with inhalation therapy.
- 3 Open the nebuliser by turning the upper part (9c) anticlockwise.
- 4 Pour the drug prescribed by the physician into the lower part (9a). Close the nebuliser by turning the upper part (9c) clockwise.
- 5 Connect the accessories as shown in the “Assembly diagram”. The upper cap (9d) provided with your “RF6 basic² nebuliser makes it possible to deliver the medication in 2 modes.


Standard mode FIG. A1

For an excellent delivery with a minimum waste of medication.

High speed mode FIG. A2

For a quicker application.

- 6 Sit comfortably, holding the nebuliser in your hands, place the mouthpiece over your mouth or use the mask. Start the unit by means of switch (5) and breathe in and out deeply.
- 7 When the treatment is finished, turn off and unplug the unit.
- 8 Clean the nebuliser and the accessories according to the instructions provided in the “Cleaning, disinfection and sterilisation” paragraph.
- 9 It may happen that after a treatment session visible humidity deposits form in the tube (8): disconnect the tube from the nebuliser and let it dry using the compressor air flow, as this prevents mould from growing in the tube.

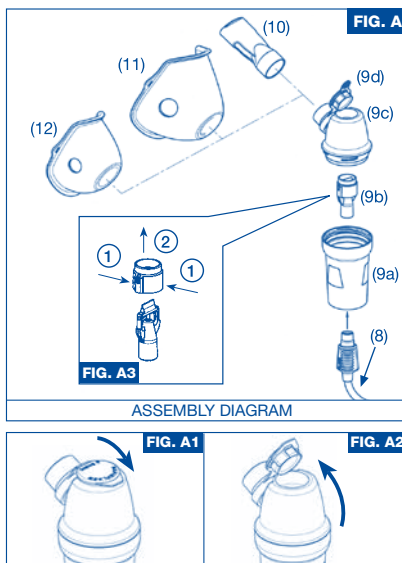2. CLEANING, DISINFECTION AND STERILISATION

Before undertaking any cleaning operation, switch off and unplug the unit.

Accessories for aerosoltherapy and Nasal wash

- **Cleaning:** Wash under running water with mild (non abrasive) washing-up liquid or in the hot cycle dishwasher, after disassembling all the parts (except for the connecting tube). Clean the nozzle (9b) after detaching it from the upper part (9c) of the nebuliser and disassembling as shown in FIG. A3.
- **Disinfection:** All accessories (except for the connecting tube) can be disinfected or sterilized choosing one of these methods: a) by soaking in a solution of water and vinegar (60% water and 40% vinegar). Finally, rinse thoroughly with drinking water; b) by boiling in water for 20 min.; this method cannot be used for masks. Use preferably distilled or demineralised water.
- **Sterilisation:** Sterilise all parts with cold chemical sterilizers (peracetic acid based, following the instructions provided by the manufacturer), or (except for tube and masks) in an autoclave at 121°C (1.2 bar) for 20 minutes, carefully following the instructions provided by the autoclave manufacturer.
- **Drying:** After washing or disinfecting the accessories, dry them very carefully with a paper tissue or, better yet, a warm air flow (e.g. from a hair dryer).
Reassemble the nebuliser nozzle (9b) and fit it in the nebuliser as indicated in the “Assembly diagram” FIG. A.

⚠ When used very frequently, we recommend that the nebuliser be replaced every 6 months (or earlier if it becomes clogged), in order to guarantee an adequate hygienic level and the maximum therapeutic effect.

Unit and exterior of the tube

To clean use only a damp cloth and an antibacterial detergent (non-abrasive and solvent free).

3. AIR FILTERING

The unit is equipped with an air filter (4) which should be replaced when dirty or when its colour changes. Do not wash or re-use the used filter. The filter must be regularly replaced to help ensure efficient compressor performance. The filter should be regularly checked. For spare filters contact your dealer or an authorised service centre.

To replace the filter:

Lift the lid and remove the filter using a screwdriver (see FIG. B). The filter is designed to remain fixed in its housing.

Use only Flaem original accessories.

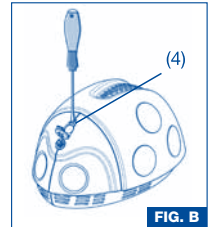
FIG. B

4. IMPORTANT SAFEGUARDS

⚠ Before using the unit for the first time, and periodically throughout its life, check the power supply cord for any damage; if you detect any damage, do not plug it in. Take the unit immediately to an authorised service centre or your dealer.

- Do not allow to be used as a toy. Close attention is necessary when used by or near children.
- Children and disabled persons should always use the unit under strict supervision of an adult who has read this manual.
- Some of the components are small enough to be swallowed by children; therefore, keep the unit out of children's reach.
- The unit should not be used in the presence of an anaesthetic mixture inflammable with air, oxygen or nitrous oxide.
- Never use adaptors for power supplies different from the one indicated on the bottom of the unit (plate details). Keep the cord away from heated surfaces.
- Use solely the provided power supply cord.
- Do not handle the plug with wet hands and do not use the unit while bathing or showering. Never submerge the unit in water; should this accidentally occur, unplug it immediately. Do not remove or touch the unit under water before pulling out the plug. Take it immediately to an authorised service centre or your dealer.
- Never wash the unit under running water or by soaking.
- The casing of the unit is not protected against the infiltration of liquids.
- Always unplug the power cord immediately after use.
- Do not obstruct the air slots placed on both sides of the unit. During use, always place the unit on a firm surface clear of obstacles. The air slots should be kept free of obstruction by blankets, pillows, hairs or similar objects. Do not place any objects inside the air slots.
- Servicing may be performed by authorised personnel only. Any unauthorised repairs will void the warranty and may pose a safety hazard for the user.
- The Manufacturer, the Retailer and the Importer are responsible for safety, reliability and performance only if: a) the unit is used in accordance with the operating instructions; b) the electrical system of the premises where the unit is used is in conformity with current laws.
- **Microbial contamination:** a personal use of the nebuliser and its accessories is recommended in order to avoid any risks of infection.
- **Interactions:** the materials used that come into contact with drugs have been tested with a wide range of medications. However, due to the variety and continuous evolution of medication, the possibility of a chemical interaction cannot be excluded. Once opened, the medication should be used as soon as possible. Avoid prolonged exposure to the nebuliser or the nasal wash.

5. TECHNICAL SPECIFICATIONS

Unit Mod.: P0406EM F400

Voltage:	230V ~ 50Hz 130VA	115V ~ 60Hz 2A Max	220V ~ 60Hz	100V ~ 50/60Hz
In conformity with:	Dir. 93/42/CEE (EN 60 601-1) 			
Approvals:				

Fuse: T 2.5A – 250V
 Max pressure: 1.5 bar approx.
 Compressor air output: 9 l/min. approx.
 Dimensions: 21x16x13 cm
 Weight: 1.73 Kg
 Sound level (at 1 m): 55 dB (A) approx
 Continuous use

RF6 BASIC² nebuliser

Medication minimum capacity: 2 ml
 Medication maximum capacity: 8 ml
 Operating pressure (with nebuliser): 0.60 bar approx.
 Max. delivery - Standard mode: 0.25 ml/min approx.
 Max. delivery - "High speed" mode: 0.32 ml/min approx.

In vitro characterization carried out at the Inamed Research GmbH & Co. KG on behalf of TÜV Rheinland Product Safety GmbH – Germany in compliance with the new European Standard for aerosol therapy appliances, Standard EN 13544-1. Further details are available upon request.

MMAD - Standard mode 3.8 µm	MMAD - "High speed" mode 4 µm
Breathable fraction < 5 µm (FPF) - Standard mode 62 %	
Breathable fraction < 5 µm (FPF) - "High speed" mode 63 %	

Operating conditions: Temperature: min 10°C; max 40°C

RH Air humidity min 10%; max 95%

Storage conditions:  Temperature: min -25°C; max 70°C

RH Air humidity min 10%; max 95%

At atmospheric pressure operating and storage conditions: **P** min 690 hPa; max 1060 hPa

 **Disposal:** In conformity with Directive 2002/96/EC, the symbol shown on the appliance to be disposed of indicates that this is considered as waste and must therefore undergo "sorted waste collection".

■ The user must therefore take (or have taken) the above waste to a pre-sorted waste collection centre set up by the local authorities, or else give it back to the dealer when purchasing a new appliance of the same type. Pre-sorted waste collection and the subsequent treatment, recovery and disposal operations favour the production of appliances made of recycled materials and limit the negative effects of any incorrect waste management on the environment and public health. The unlawful disposal of the product by the user could result in administrative fines as provided by the laws transposing Directive 2002/96/EC of the European member state in which the product is disposed of.

SYMBOLS

 Class II unit

 Type B unit

 Important safeguards

 Important: check the operating instructions

 Alternating current

 Switch OFF

 Switch ON

 Do not use the unit while bathing or showering

 Latex Free

Electromagnetic compatibility:

The unit is designed to meet the current requirements to meet the current requirements for electromagnetic compatibility. However if you suspect that the appliance is interfering with the operation of your TV, radio or other electrical appliance, try relocating the unit until the interference ceases. Alternatively try plugging the appliance into a different power socket.

The Manufacturer reserves the right to make technical and functional changes to this product without prior notice.

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